



CLINIPET+

single channel fixed volume pipettes




SCIENCE AT YOUR FINGERTIPS



CLINIPET[†] single channel fixed volume pipettes

Throughout years the major advantage of CLINIPET⁺ fixed volume pipettes were their outstanding repeatability, ergonomics, combined with a very competitive price. Unlike most of fixed volume pipettes, CLINIPET⁺ is not based on the variable volume pipette handle. The CLINIPET⁺ handle and mechanism was designed especially for the fixed volume pipette range, consequently not possessing the components affecting pipette weight or affordability. The pipette features reduced pipetting and blowout forces, highly durable shafts and a fully autoclavable, UV resistant body. CLINIPET⁺ is an ideal pipette for a number of clinical and educational applications.

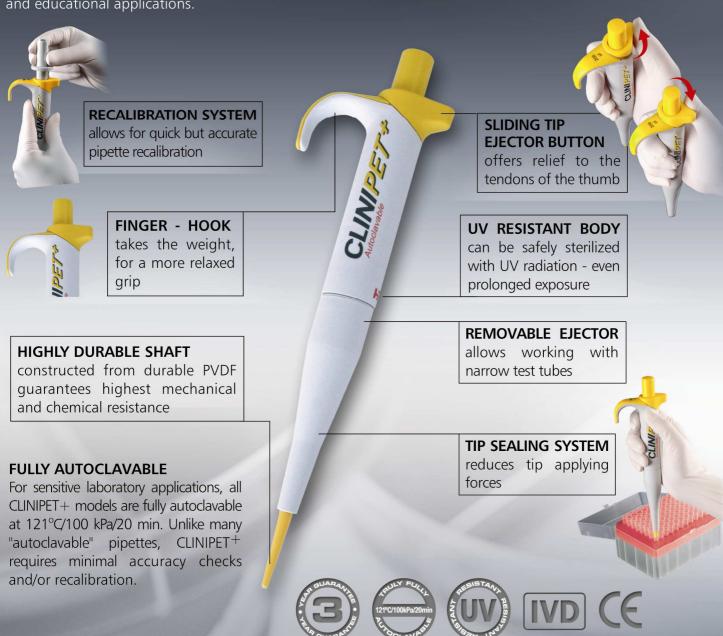

Revolutionary tip sealing system

CLINIPET⁺ pipette is equipped with a unique, patented tip sealing system. The shaft of CLINIPET⁺ retrieves when it is pressed against the pipette tip. Thanks to this feature the force required for tip sealing is reduced to minimum. Safely sealed, with the lowest possible force, the tip can be ejected with virtually no effort.

CLINIPET⁺ pipette series ordering information


Model	Cat. no.	Colour code	Volume (μl)	A(%)	P(%)	Fit to tips (μl)
CP2	2202		2	± 3.0	± 1.5	200
CP5	2204		5	± 2.0	± 1.0	
CP7	2206		7	± 1.5	± 0.7	
CP10	2208	Yellow	10	± 1.2	± 0.5	
CP15	2210		15	± 1.0	± 0.5	
CP20	2212		20	± 0.9	± 0.4	
CP25	2214		25	± 0.9	± 0.4	
CP30	2216		30	± 0.9	± 0.3	
CP40	2218		40	± 0.9	± 0.3	
CP44.7	2220		44.7	± 0.9	± 0.3	
CP50	2222		50	± 0.9	± 0.3	
CP60	2224		60	± 0.9	± 0.3	
CP70	2226		70	± 0.9	± 0.3	
CP75	2228		75	± 0.9	± 0.3	
CP80	2230		80	± 0.9	± 0.3	
CP90	2232		90	± 0.8	± 0.3	
CP100	2234		100	± 0.8	± 0.3	
CP120	2236		120	± 0.8	± 0.3	
CP150	2238		150	± 0.7	± 0.3	
CP200	2240		200	± 0.6	± 0.2	
CP200A	2242		200	± 0.6	± 0.2	1000
CP220	2244		220	± 0.6	± 0.2	
CP250	2246		250	± 0.6	± 0.2	
CP300	2248		300	± 0.6	± 0.2	
CP400	2250		400	± 0.6	± 0.2	
CP450	2252		450	± 0.6	± 0.2	
CP500	2254	Blue	500	± 0.6	± 0.2	
CP600	2256		600	± 0.6	± 0.2	
CP700	2258		700	± 0.6	± 0.2	
CP750	2260		750	± 0.6	± 0.2	
CP800	2262		800	± 0.6	± 0.2	
CP900	2264		900	± 0.6	± 0.2	
CP1000	2266		1000	± 0.6	± 0.2	





Pipette stands

In order to ensure convenient storage of pipettes, HTL offers a wide range of single and multiple stands.

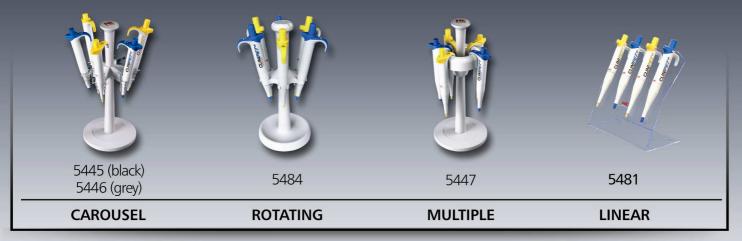

HTL pipette tips

HTL Pipette Tips are designed to fit a wide variety of single and multichannel pipettes. Flexible walls and the diameter of tips collars ensure a secure fit with less force required to load and eject. All pre-sterilized tips are exposed to e-beam radiation for complete security without degrading the plastic or performance. Tips are manufactured in one of the most technically advanced Axygen plastic disposable production facility in clean room conditions. HTL Pipette Tips are certified to be free from bio active contaminants as RNase, DNase or pyrogens. HTL tips are offered in standard, sterile and filtered versions in all popular packaging formats:



Distributed by: